


## Mammoth DETECTOR COVID-19 Test summary

Notebook: tsdobbs's notebook

Created: 3/26/2020 11:38 PM

Author: Tim Dobbs

Updated: 3/27/2020 3:51 AM

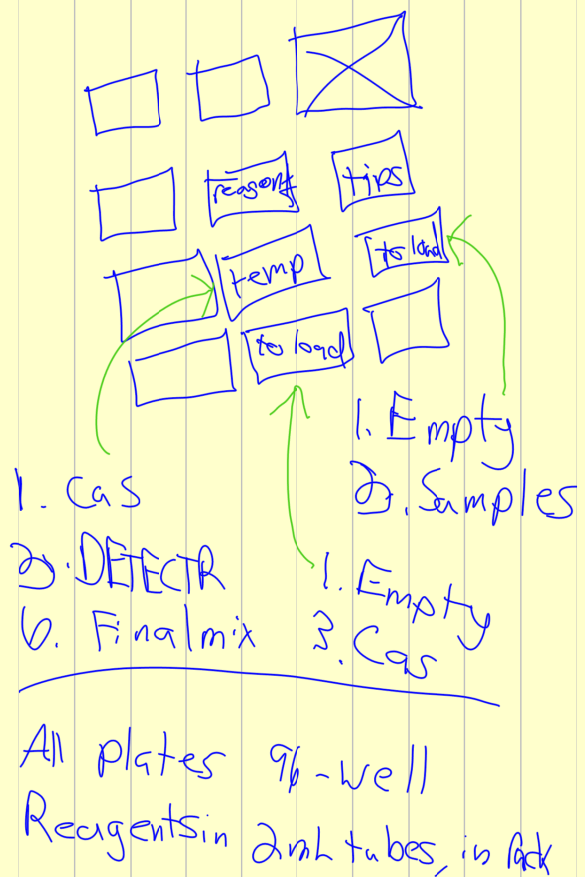

- 1) Prep Lbcas12on automatically, in batches.  
chill
- 2) user moves to ice then sets up new lab ware
- 3) In still-cold block Prep 3 aliquots for each sample
- 4) Incubate @  $62^{\circ}\text{C}$  for 20 min  
Cool to  $37^{\circ}\text{C}$
- 5) User moves pre-prep Lbcas12on to deck
- 6) Mix. Notify user after 10 min